



# Bionic Board Autonom Concept



[www.beamyourlife.com](http://www.beamyourlife.com)  
Email: [info@beamyourlife.com](mailto:info@beamyourlife.com)  
Tel: +49 (0)8024 82 63  
Made in Germany

**Manual**  
**Model 2024**



## Guarantee

The guarantee is valid for 2 years after the purchase of the product for private use. The general terms and conditions apply to the warranty.

Please keep the invoice, the manual and the warranty carefully.

## Consumer rights

The place of performance is the registered office of the manufacturer or importer as stated in the manual. If the damage was caused by an operating error or incorrect operation by the customer or by natural disasters, service can only be provided for an additional fee.

Especially in the following cases, the expenses incurred will be invoiced:

- For damage caused by incorrect operation
- For damage caused by using the product in conditions not recommended
- For damage to wearing parts
- For damage caused by a blow or by the product falling on a hard object
- In case of damage caused by connection to the wrong mains connection or by overvoltage
- For damage caused by not following the manual
- Store and use the device only in closed and dry rooms
- Do not use it near liquids or flammable gases
- Do not attempt to repair, open or disassemble the unit yourself

Thank you for choosing the Bionic Board.

Please read this manual carefully so that you can use the advantages of the product safely and effectively.

Please keep this manual in a safe place so that you can refer to it later if necessary.

## Technical data

Voltage: AC 220 - 240 Volt, 50 Hz  
Power: approx. 200 Watt  
Weight: approx. 6,5 kg  
L x B x H: 49 x 31 x 11 cm

## Table of contents

### Contraindication

### Safety Instructions

### Operation

### Guarantee

### Consumer rights

#### **Note!**

In the rare case that the product does not function properly after switching on, please remove the power plug from the socket for about 10 seconds and then switch the product on again.

Do not use the appliance for more than 15 minutes in continuous operation to avoid overheating!

Let the product cool down sufficiently before using it again.

Ensure adequate ventilation during operation and do not cover the product.

## Contraindication

The following groups of people should not use the device or consult a doctor before using the product:

- Pregnant women
- Persons with (infectious) diseases of all kinds
- Persons with electronic pacemakers of all kinds

This product should not be used on burns.

## Safety Instructions

To avoid the danger of electric shocks:

Always remove the plug of the power supply unit from the socket immediately after use and before cleaning.

Do not use the appliance in temperatures above 40 degrees Celsius.

Keep the product dry. Do not use the product in damp or wet environments.


If the power cord is damaged, it must be replaced immediately by the manufacturer, a service technician or similar qualified person to avoid the risk of electric shock.

The product was developed for domestic use and not for commercial use.

To reduce the risk of fire, electric shock or injury to persons, connect the product only to a suitable wall outlet and the supplied AC adapter.

If the product is moved from a cold environment to a warm environment, wait at least 1 hour before using the product, so that possible condensation in the electronic parts does not cause a short circuit.

## Operation

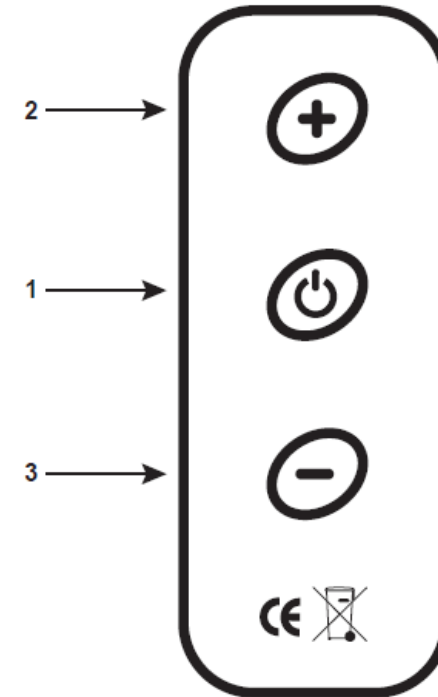

- 1 ON / OFF switch
- 2 "+" intensity higher
- 3 "-" Intensity lower

Press button 1 (ON / OFF) to start. With buttons 2 and 3 you can adjust the intensity in 5 steps.

The duration of the automatic program is about 10 minutes.

For an optimal mobility effect, please also move your feet horizontally and vertically during the usage of the product.